

OXLIFE INDEPENDENCE®

User Manual



PORTABLE OXYGEN CONCENTRATOR

Oxlife Independence® 301-0001V and 301-0001CE

USER: Read this manual before operating this device. Save this manual for future reference.

HEALTH-CARE PROVIDER: This manual must be provided to the end user.

Table of Contents

- Section 1: Introduction 5**
 - Intended Purpose 6**
 - Intended Users 6**
 - Intended Patient Population 6**
 - Indications 6**
 - Contraindications..... 6**
 - Device and Accessory Symbol Reference 7**
 - User Interface Symbols..... 8**
 - Specifications..... 9**
 - DNA Technology® 11**
- Section 2: Safety Guidelines 12**
 - Device Safety Guidelines..... 12**
 - Patient Safety Guidelines 14**
 - Battery and Power Supply Safety Guidelines..... 16**
- Section 3: Product Description 17**
 - Device and Accessories 17**
 - Feature Identification 18**
 - Understanding the Control Panel 19**
- Section 4: Operating Instructions 20**
 - Before Operating 20**
 - Locating Your Device..... 20**
 - Device Settings..... 21**
 - Installing and Removing the Battery 23**
 - Typical Battery Operation Times..... 23**
 - Battery Time Management..... 24**
 - AC Power Supply 25**
 - DC Power Cord 26**
 - DC Power Cord User Information Guide 27**
 - DC Power Cord Troubleshooting Guide 27**

Handle Operation	28
Accessory Bag.....	28
Cannula Use.....	29
Accessing the Provider Screen	30
Changing Languages.....	30
Section 5: User Alerts and Alarms.....	32
Alarm and Alert Screens.....	32
Alarm and Alert Screen Descriptions	33
Alarm System Test	37
Section 6: Maintenance and Cleaning.....	38
Maintenance	38
Pre-Use Functional Check	38
Device Care and Cleaning	39
Device Disposal.....	40
Battery Disposal	40
Section 7: Traveling	41
Entering and Exiting Airplane Mode	41
Section 8: Standards Compliance	42
Classification	42

Section 1: Introduction

PLEASE READ THIS OPERATION MANUAL CAREFULLY BEFORE USING THIS DEVICE. BE AWARE OF ALL WARNINGS AND SAFETY INFORMATION. ONLY USE ACCESSORIES APPROVED BY O2 CONCEPTS AND REFERENCED WITHIN THIS MANUAL.

IF YOU DO NOT FULLY UNDERSTAND ALL THE WARNINGS, SAFETY PRECAUTIONS, AND OPERATING INSTRUCTIONS CONTACT YOUR AUTHORIZED DEALER OR HEALTHCARE PROVIDER FOR TECHNICAL SUPPORT.

Please contact your authorized dealer or healthcare provider if your Oxlife Independence[®] requires service.

Please call O2 Concepts[®] Technical Support to report any unexpected operation or events associated with the device at 1-877-867-4008 EXT 367 or contact the EU Rep or other Authorized Representative listed on the rear cover of this user manual, or your country's relevant government agency, such as the competent authority.

Additional information about O2 Concepts[®], the Oxlife Independence[®], and the device limited warranty, can be found on our website at www.o2-concepts.com.

CAUTION: U.S. FEDERAL LAW RESTRICTS THESE DEVICES TO SALE BY ORDER OF A PHYSICIAN.

Intended Purpose

The Oxlife Independence® is intended to administer supplemental oxygen delivered nasally to adult and pediatric patients. It is not intended for life support.

Intended Users

The Oxlife Independence® is intended to be used by patients in a manner according to the intended purpose and operating guidelines. These devices are designed to be operated by patients who have received appropriate supervision and training, as needed, from healthcare professionals to ensure safe and effective use.

Intended Patient Population

The patient population includes pediatric and adult patients requiring supplemental oxygen as prescribed by a healthcare provider in the treatment of disorders aligning with the indications and contraindications of the device.

Indications

The Oxlife Independence® is indicated on a prescription basis for the administration of supplemental oxygen used in the treatment of hypoxemia and other related conditions in patients with Chronic Obstructive Pulmonary Disease (COPD), Interstitial Lung Disease (ILD), Congestive Heart Failure (CHF), or other disease conditions where the administration of supplemental oxygen may be appropriate. It is not intended for life support. It is a transportable device capable of continuous use in the home, institutional, and vehicle/mobile environments.

Contraindications

- ⚠ Under certain circumstances, the use of non-prescription oxygen therapy can be hazardous. This device should only be used when prescribed by a healthcare provider.
- ⚠ Not for use in the presence of aerosol sprays or flammable anesthetics.
- ⚠ Additional monitoring may be required, as appropriate, for use by geriatric, pediatric, or any other patient unable to hear, see or interpret alarms, or communicate discomfort.
- ⚠ The Oxlife Independence® is not appropriate for any patient who would experience adverse health consequences as a result of a temporary interruption in oxygen therapy.
- ⚠ The availability of an alternate source of supplemental oxygen is recommended in the case of a power interruption or a mechanical failure of the device. Consult your healthcare provider for their recommendation.
- ⚠ It is the patient's responsibility to arrange for an alternate oxygen supply when traveling. O2 Concepts assumes no liability for persons choosing not to adhere to manufacturer recommendations.
- ⚠ Use of the device by patients with tracheotomy, under the recommendation of a qualified healthcare provider, will require the use of appropriate clinical apparatus or attachments.

Device and Accessory Symbol Reference

ICON	MEANING	ICON	MEANING
DEVICE SYMBOLS			
	Caution represents the possibility of damage to the equipment.		Warning represents the possibility of harm to the operator or patient
 	No smoking while using or near device		European Declaration of Conformity
 	No open flames		Medical Device
 	Refer to instructions for use	 	Use no grease or oils
 	Do not disassemble		Drip proof equipment
	Not to be serviced by users		Keep dry in transport and storage
	Type BF equipment		Manufacturer
	Prescription Use Only		Class II equipment
ACCESSORY SYMBOLS			
	Recycle battery		Indoor Use Only*
	Do not incinerate battery		Do not dispose of in household waste
	AC Power (Alternating current)		DC Power (Direct Current)

Table 1 *Specific to the AC power supply and DC power cord, not the Oxlife Independence device.

User Interface Symbols

ICON	MEANING	ICON	MEANING
	Power button		Breath detection
	Mode button		Battery life indication
	Continuous Mode		AC/DC power
	Pulse Mode		AC/DC charge indication
	Increase flow setting Scroll Up		Contains DNA Technology®*
	Decrease flow setting Scroll Down		

Table 2 *The DNA Technology label is included on enabled devices.

Specifications

DIMENSIONS		
Device Dimensions with Handle/Wheels	H: 20.29 in W: 10.85 in D: 9.45 in	H: 51.5 cm W: 27.55 cm D: 24 cm
WEIGHTS		
Device Weight	16.7 lbs. (7.57 kg)	
Cart Kit (Wheels & Pull Handle)	2.17 lbs. (0.98 kg)	
Battery	1.4 lbs. (0.63 kg)	
AC Power Supply	1.9 lbs. (0.86 kg)	
DC Power Cord	.37 lbs. (0.16 kg)	
MODES OF OPERATION		
Continuous Flow	Six settings; 0.5 to 3 LPM in 0.5 LPM increments; Maximum deviation in flow rate is $\pm 10\%$. Max flow 3.0 LPM. Flow maintained with outlet pressure ranging from 0 to 5.0 psig (0-34kPa)	
Pulse Dose	Twelve settings; 0.5 to 6.0 (8mL-96mL) in 8ml increments; Maximum deviation of volume per breath is $\pm 15\%$ across all environmental conditions.	
Breath Sensitivity	The device will deliver bolus based on a trigger pressure range of -0.08 to -0.35 cmH2O.	
Battery Specification	14.4V Lithium Ion Battery	
External Power Supply Power Input	AC power: 100 - 240 VAC, 50/60 Hz, 2.0A DC power: 12-15 VDC; Recommended: 10A outlet at 12V	
Operating Altitude	0-13,123ft or (0-4000m)	
Oxygen Purity	91% $\pm 4\%$ (87 - 95%) from 0.5 to 3.0 LPM; Measured purity values are within $\pm 2\%$ of actual values; There is no variation in purity within the operating altitude. This applies to the full range of environmental operating conditions.	
Operating Temperature	50°F (10°C) to 104°F (40°C)	
Operating Humidity	10% - 95% @ 82.4°F (28°C)	
Operating Atmospheric Pressure	101kPa to 63kPa	
Operating Environment	Free of smoke, pollutants, and fumes.	
Transport/ Storage Temperature	-4°F (-20°C) to 140°F (60°C)	
Transport/ Storage Humidity	0 -95% non-condensing	
Operating Time	24 hours per day when connected to an external AC or DC power source. This is a continuous operation device.	
Cannula Specification	DO NOT use cannula tubing longer than 7ft (2m) when using Pulse Flow Mode Do NOT use cannula tubing longer than 50ft (15m) when using Continuous Flow Mode.	

Sound Level	The typical volume at Pulse Setting 2 will be ~38.6 dBA. The peak volume the device will reach at maximum settings will be ~46.6 dBA.
Alarm Sound Level	All alarms triggered by the device will be 85dBA or higher at 10cm from the device.
Maximum Outlet Pressure	Device will maintain flow up to 5psi back pressure. Maximum back pressure is 7.0psi.
Service Life	Device: 5yrs, Accessories: 1yr

DNA Technology®

Some Oxlife Independence® devices contain DNA Technology® (Dynamic Network Analysis), which is available to allow your health care provider to better serve your needs. This technology is intended to assist your provider in determining the device is operating properly.



This device is FAA approved for use aboard passenger aircraft (FAA Advisory Circular 91.21-1B), including radio frequency emission limits of (RTCA) Document (DO) 160. section 21. Category M. Device contains FCC Id XPY2AGQN4NNN.

The cellular connection must be deactivated prior to flight on any commercial aircraft.

Instructions for deactivating the cellular connection are listed on page [41](#) within Section 7: Entering and Exiting Airplane Mode.

There is no known interference posed by medical equipment during specific investigations or treatments.

There are no known devices that will cause interference issues.

Applied Parts

The applied parts, or components that a patient will come in contact with during normal operation of the device, are listed below:

- Pull and Carry Handles
- Control Panel
- Oxygen Outlet Port
- External Power Input

Section 2: Safety Guidelines

Device Safety Guidelines

- ❗ AVOID EXPOSURE TO OPEN FLAMES OR CREATION OF ANY SPARK NEAR YOUR OXLife INDEPENDENCE®. THIS INCLUDES SPARKS FROM STATIC ELECTRICITY CREATED BY ANY TYPE OF FRICTION. PROTECT ELECTRICAL POWER CORDS FROM SHARP EDGES TO AVOID ELECTRICAL SHOCK AND SERIOUS PHYSICAL INJURY.
- ❗ Locate the Oxlife Independence® in a well-ventilated area to allow for adequate air intake. Avoid the intake of airborne pollutants, smoke, or fumes. Wind or strong drafts can adversely affect accurate delivery of oxygen therapy.
- ❗ Use only approved accessories as specified in this User Manual to ensure proper function and avoid the risk of fire and burns. Use of non-approved accessories may cause serious damage to the device and will void the warranty.
- ❗ Locate oxygen tubing and power cords away from hot surfaces and in a manner to prevent tripping hazards.
- ❗ **DO NOT** modify the device. Any modifications made to the device may impair performance or damage the equipment and will void the warranty.
- ❗ **DO NOT** operate the device in an enclosed space.
- ❗ **DO NOT** cover the device or block the air inlet or the exhaust ports as it may impair device function.
- ❗ **DO NOT** drop or insert objects or liquid into any opening.
- ❗ **DO NOT** leave your Oxlife Independence® or batteries in your vehicle or trunk. Extreme heat or cold may damage your device and/or batteries.
- ❗ When shipping the Oxlife Independence® and accessories, the batteries must be wrapped separately to avoid accidental battery contact. **DO NOT** ship the Oxlife Independence® with the batteries installed.
- ❗ Keep the device away from potential household pests to avoid infestation that will impair device performance and will void the warranty.
- ❗ This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:
 1. This device may not cause harmful interference, and
 2. This device must accept any interference received, including interference that may cause undesired operation.
- ❗ Portable RF communications equipment (including peripherals such as antenna, cables, and external antennas) should be used no closer than 12 inches (30cm) to any part of the Oxlife Independence® including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

- ① Use of accessories, transducers, or cables other than those specified by O2 Concepts could result in increased electromagnetic emissions or decreased electromagnetic immunity of the device and result in improper operation.
- ① Use of the device adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, the device as well as other stacked equipment should be observed to verify that they are operating normally.
- ① The device is suitable to use in a CISPR 11 Group 1 Class B environment for radiated emissions.

Patient Safety Guidelines

- ⚠ PRESCRIBED MODE AND FLOW SETTINGS SHOULD ONLY BE ADJUSTED UNDER THE ADVICE OF A HEALTHCARE PROVIDER.
- ⚠ If you feel discomfort or are experiencing a medical emergency while undergoing oxygen therapy, seek medical assistance immediately to avoid harm.
- ⚠ To ensure you are receiving the therapeutic amount of oxygen according to your medical condition, the Oxlife Independence must be used only after one or more settings have been individually determined or prescribed by a healthcare provider for you at your specific activity levels. The Oxlife Independence® must be used with the specific combination of parts and accessories that were used while your settings were determined.
- ⚠ KEEP ALL MATCHES, LIT CIGARETTES, CANDLES, OR OTHER SOURCES OF IGNITION AT LEAST 10 FEET (3M) FROM THE DEVICE. THIS DEVICE PRODUCES ENRICHED OXYGEN GAS WHICH ACCELERATES COMBUSTION.
- ⚠ **DO NOT SMOKE WHILE USING THIS DEVICE. SMOKING DURING OXYGEN THERAPY IS DANGEROUS** and is likely to result in facial burns or death. **DO NOT** allow smoking within the same room where the oxygen concentrator or any oxygen carrying accessories are located. If you intend to smoke, you must always turn the oxygen concentrator off, remove the cannula and leave the room where the oxygen concentrator is located. If you are unable to leave the room, you must turn the oxygen concentrator off and remove the cannula before smoking.
- ⚠ Use only water-based lotions or salves that are oxygen compatible before and during oxygen therapy. **DO NOT** use oil, grease, or petroleum-based products on or near the device to prevent accidental ignition.
- ⚠ **DO NOT** lubricate fittings, connections, tubing or other accessories to avoid the risk of fire and burns.
- ⚠ **DO NOT** leave the nasal cannula or mask on bed coverings or chair cushions. If the oxygen concentrator is turned on, but not in use, the oxygen will make the materials flammable. Turn the oxygen concentrator off when not in use to prevent oxygen enrichment.
- ⚠ **DO NOT** operate the device on wet surfaces or in standing water, and do not submerge or expose to water or precipitation to avoid serious injury or damage to the device. Do not contact the device when wet to avoid serious injury or the chance of shock. If the Oxlife Independence has been dropped, damaged, or exposed to water please contact your authorized dealer or healthcare provider for inspection and possible service of the device.
- ⚠ **DO NOT** use the device with an extension cord.
- ⚠ **DO NOT** use the device with a damaged power cord or plug to avoid injury.
- ⚠ Outdoor use must be conducted using battery power.
- ⚠ Replace cannula as recommended by the cannula manufacturer. Cannulas may become contaminated with body fluids or expired gases during use.
- ⚠ Monitor children and pets when in proximity to the device. Keep cannula tubing

away from children and pets to avoid danger of choking or strangulation.

- ⚠ Use of this device outside of the stated altitude, temperature or humidity ranges is expected to adversely affect the flowrate and the percentage of oxygen and, consequently, the quality of oxygen therapy.

Battery and Power Supply Safety Guidelines

- ⚠ Use of a non-approved battery or power supply may cause serious damage to the device and will void the warranty.
- ⚠ Remove the battery from the device if it will not be used for an extended period to enhance battery life. Store the battery in a cool dry place when not in use.
- ⚠ Locate the external power supply in the open air to prevent overheating when in use.
- ⚠ U.S. Department of Transportation (DOT) and United Nations (UN) regulations require that batteries be removed from the device when checked as luggage on international flights.
- ⚠ The battery may explode if exposed to or disposed of in a fire.
- ⚠ Keep the battery away from children to prevent injury.
- ⚠ **DO NOT** drive over, drag or place objects on the power cords to avoid damage.
- ⚠ Use of a frayed or damaged AC power supply or DC power cord may compromise basic safety.
- ⚠ **DO NOT** attempt to disassemble the battery or power supply. Doing so will void the warranty and may cause personal injury.
- ⚠ **DO NOT** use the battery or power supply for anything other than its intended purpose. Doing so may damage the device or cause personal injury.
- ⚠ **DO NOT** drop the battery or expose it to mechanical shock. If the battery is damaged, discontinue use and dispose of it properly to avoid personal injury.
- ⚠ **DO NOT** expose the battery to water or other liquids. If the battery gets wet, discontinue use, and dispose of it properly to avoid personal injury.
- ⚠ **DO NOT** expose the battery to excessive heat or cold outside of specifications as it may affect performance.
- ⚠ Use caution when handling the DC plug adapter. This plug may get hot with use. Ensure DC plug socket is clean of debris which may cause overheating.

Section 3: Product Description

Device and Accessories

The Oxlife Independence® is a portable oxygen concentrator used on a prescriptive basis that enables patients requiring supplemental oxygen to be treated in a home, institutional, or vehicle/mobile environment.

The Oxlife Independence® delivers 87%-95% pure oxygen to a patient through a standard single lumen nasal cannula in both continuous and pulse mode.

Pulse mode detects a patient's breath and then delivers a bolus of oxygen during inhalation. In continuous mode the device delivers a continuous flow of oxygen. Refer to the Specifications and other sections of this manual for details on specific settings and flow rates.

The Oxlife Independence® has standard power options which include a 100- 240V AC power supply, a 10-15V DC power cord, and a rechargeable lithium-ion battery.

Certain models may come equipped with O2 Concepts® DNA Technology® to provide equipment performance, usage data and GPS/GNSS location to O2 Concepts® via cellular connectivity.

- ❗ The use of certain humidifiers and administration accessories not specified for use with this oxygen concentrator may impair performance. Your healthcare provider or authorized dealer will recommend the proper accessories for your use.
- ❗ Your healthcare provider or authorized dealer will recommend the proper cannula for your use. Ensure that the selected cannula is compliant with ISO 80601-2-69 Medical Electrical Equipment, Particular requirement for basic safety and essential performance of oxygen concentrator equipment.
- ❗ The configuration of the equipment and accessories must be determined for each individual patient.


This manual applies to the following accessories:

- › Battery
- › AC Power Supply
- › DC Power Cord
- › Accessory Bag
- › Firesafe™ Cannula Valve
- › Cannula (not included)

This manual does not apply to the following accessory items sold separately:


- › Battery Shell/Blank
- › Desktop Battery Charger, Single Bay
- › Desktop Battery Charger, Dual Bay

Feature Identification

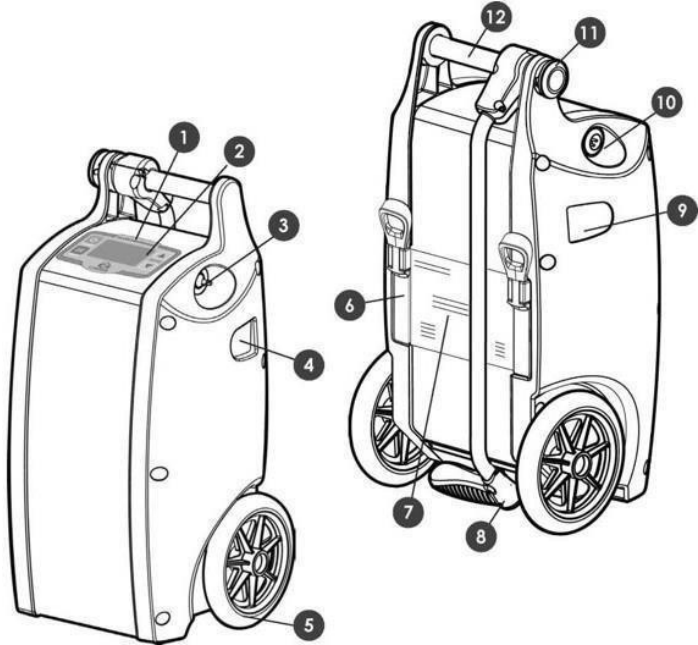

- 1. Control Panel
- 2. LCD Display
- 3. Oxygen Outlet Port
- 4. Exhaust Vent
- 5. Wheels
- 6. Rechargeable Battery Slots
- 7. FAA Approval Identification
- 8. Pull Handle
- 9. Air Intake Vent Filter
- 10. External Power Input
- 11. Handle Lock/Release Button
- 12. Carry Handle

Understanding the Control Panel

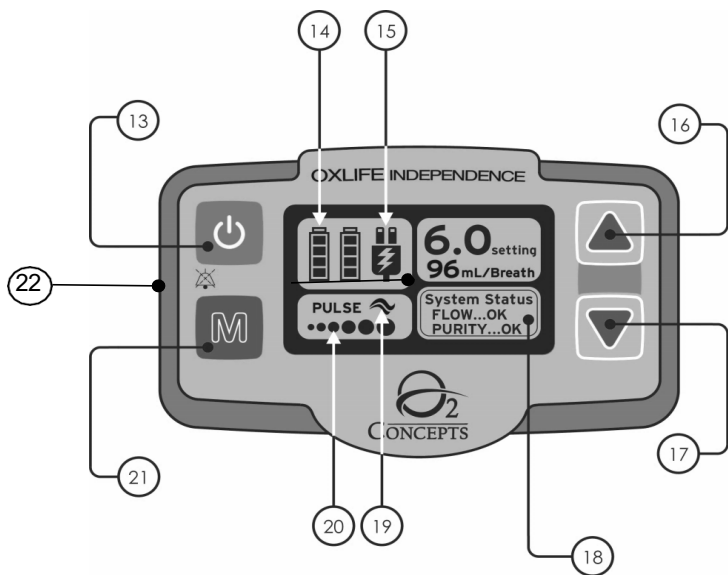

- 13. Power Button
- 14. Battery Status
- 15. External Power Source/Battery Charging Indicator
- 16. Increase/Scroll Up Button
- 17. Decrease/Scroll Down Button
- 18. System Status
- 19. Breath Detect (Pulse Mode Only)
- 20. Mode Setting
- 21. Mode Button
- 22. Indefinite Acknowledge (Alarm Mute)

Section 4: Operating Instructions

Before Operating

Your authorized dealer has demonstrated the proper operation of your Oxlife Independence®. This manual provides product information and operating instructions and should be saved for future reference to help you safely operate your device. If you have any further questions, contact your authorized dealer or healthcare provider.

This equipment must be placed into service in accordance with the information provided in the accompanying documents.

DO NOT operate the Oxlife Independence® without first reading the Safety Guidelines included in **SECTION 2** of this manual.

Please follow all operating instructions.

Before each use, ensure that filters are clean and in place and that cannula and power cord connections are intact.

If you are relocating your Oxlife Independence® from an extreme environment, allow the device to return to the specified operating temperature and humidity ranges before use.

Operating your device outside of specified ranges may damage your device, impact device performance and may void your warranty. See Operating Temperature and Humidity Ranges listed in the Specification Table.

Fully charge batteries before first use.

The essential performance of the Oxlife Independence® is to produce $91\% \pm 4\%$ oxygen at all continuous and pulse settings.

Locating Your Device

Proper placement and positioning of the device is critical to the effectiveness of oxygen therapy.

Place the Oxlife Independence® in a well-ventilated area free of smoke, fumes, pollutants, and away from direct sunlight. Avoid high humidity environments. Ensure that air intake and exhaust ports are not obstructed, and air intake vent filters are in place.

The Oxlife Independence® **MUST** be located so that alarms can be heard. Position the oxygen supply tubing and power cords in a manner that prevents kinking, air flow obstructions, and tripping hazards.

DO NOT locate the Oxlife Independence® near any flammable materials or cleaning product or in the direct path of any heat source such as a stove, heat register or a car heater. Keep the Oxlife Independence® at least 10 feet (3m) from hot sparking objects or open flame.


Device Settings

Turn the device on by pressing the Power Button .

Press the Mode Button  to select Continuous or Pulse Flow Mode.

Press the Increase  or Decrease  buttons to select the correct flow rate.

Turn the device off by pressing and holding the Power Button .

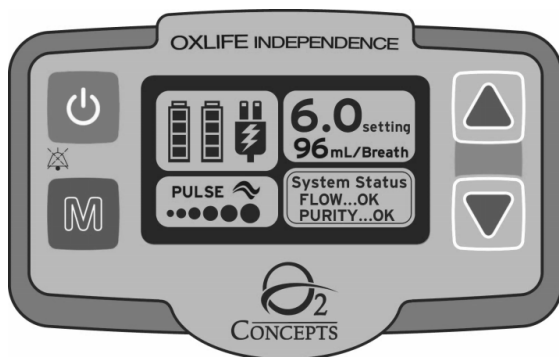

The device may take a minimum of seven (7) minutes to warm up and reach the desired performance. System Status will show Flow and Purity are OK when warm-up is complete.

THE PROPER FLOW MODE AND FLOW RATE ARE PRESCRIBED BY YOUR HEALTHCARE PROVIDER. DO NOT CHANGE THESE SETTINGS WITHOUT CONSULTING YOUR HEALTHCARE PROVIDER.

PULSE FLOW MODE SHOULD ONLY BE USED UNDER THE DIRECTION OF YOUR HEALTHCARE PROVIDER. APPROPRIATE AND SAFE PULSE MODE SETTINGS MUST ACCOMMODATE THE INDIVIDUAL PATIENT'S LIFESTYLE INCLUDING REST, TRAVEL AND EXERCISE.

DEVICE OXYGEN DELIVERY SETTINGS SHOULD BE PERIODICALLY REASSESSED FOR THE EFFECTIVENESS OF THE OXYGEN THERAPY.

THE PATIENT IS THE INTENDED OPERATOR.

The Oxlife Independence® is provided on a prescription basis for the administration of supplemental oxygen. Supplemental oxygen may reduce breathlessness and improve your ability to maintain an active lifestyle.

Consult your healthcare provider if you are feeling unwell, which may indicate too much or too little oxygen. These side effects are not immediate or life threatening. To monitor oxygen levels, your healthcare provider may recommend the use of a pulse oximeter.



CONTINUOUS FLOW MODE

In Continuous Mode the device will deliver a continuous flow of oxygen through the nasal cannula. The oxygen is measured in Liters per Minute or LPM.



PULSE FLOW MODE


In Pulse Flow Mode the device will detect your breath and deliver a measured pulse of oxygen or bolus. The breath detect icon will flash on the control panel to indicate breathing is detected.

Pulse Mode Bolus Volumes

Breaths per Minute (BPM)	Setting 1	Setting 2	Setting 3	Setting 4	Setting 5	Setting 6	Setting 7	Setting 8	Setting 9	Setting 10	Setting 11	Setting 12
	0.5	1.0	1.5	2.0	2.5	3.0	3.5	4.0	4.5	5.0	5.5	6.0
15	8.0	16.0	24.0	32.0	40.0	48.0	56.0	64.0	72.0	80.0	88.0	96.0
20	8.0	16.0	24.0	32.0	40.0	48.0	56.0	64.0	72.0	80.0	88.0	96.0
25	6.4	12.8	19.2	25.6	32.0	38.4	44.8	51.2	57.6	64.0	70.4	76.8
30	5.3	10.7	16.0	21.3	26.7	32.0	37.3	42.7	48.0	53.3	58.7	64.0
35	4.6	9.1	13.7	18.3	22.9	27.4	32.0	36.6	41.1	45.7	50.3	54.9
40	4.0	8.0	12.0	16.0	20.0	24.0	28.0	32.0	36.0	40.0	44.0	48.0

- ⚠ Shallow respiratory effort of the patient may not trigger a bolus.
- ⚠ The pulse flow settings of the Oxlife Independence® do not correspond with continuous flow oxygen.
- ⚠ The settings of other models or brands of oxygen therapy equipment do not correspond with the settings of the Oxlife Independence®.

Installing and Removing the Battery

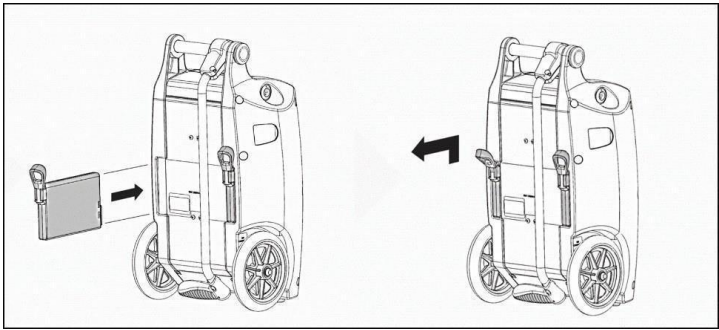

INSTALL	REMOVE
SLIDE INTO BACK OF DEVICE AND PRESS TO ENGAGE; BATTERY WILL DROP DOWN SLIGHTLY AND CLICK INTO POSITION	PULL TAB OUT SLIGHTLY TO RELEASE IT FROM THE BAY; PULL UP AND OUT ON TAB AND THEN SLIDE STRAIGHT OUT

Typical Battery Operation Times

FLOW RATE SETTING	CONTINUOUS FLOW (Hours of Battery Life)		PULSE FLOW (Hours of Battery Life)	
	1 Battery	2 Batteries	1 Battery	2 Batteries
0.5	2.4	4.75	2.5	5
1.0	1.9	3.75	2.25	4.5
1.5	1.5	3	2.1	4.25
2.0	1.25	2.25	1.9	3.75
2.5	1	2	1.75	3.5
3.0	0.75	1.5	1.5	3
3.5	Times shown are based upon 20 breaths a minute		1.4	2.75
4.0			1.25	2.5
4.5			1	2
5.0			1	2
5.5			0.9	1.75
6.0			0.9	1.75

Battery operating times are based on a new, fully charged battery and are +/- 10%.

Battery operating time will degrade with battery age, number of charge cycles and operating environment.

Battery Time Management

The Oxlife Independence® comes equipped with a rechargeable lithium-ion battery that is not user serviceable.


The batteries utilized will degrade over time from standard use. The rate of such degradation is dependent upon the frequency of operation of the device under battery power as well as frequency and length of associated battery charging cycles.

The manufacturer recommends that the user regularly monitors battery performance by viewing the battery charge level indicator on the battery and the device display. Batteries should be replaced at least every two years.

With the combined use of the battery, AC power supply, and DC power cord you may extend your time away from home. Using the AC and DC power cords whenever possible will improve battery life. The battery will charge in the device, when plugged into an external power source.

Battery Charge Status is displayed on the control panel. Each bar represents approximately 25% of the total battery charge. When the battery is fully charged (over 90%) the battery icon will appear solid.

The Oxlife Independence® battery includes a battery charge status indicator located on the front of the battery. Simply press the "PUSH" button on the battery to display the remaining battery life in 25% increments.

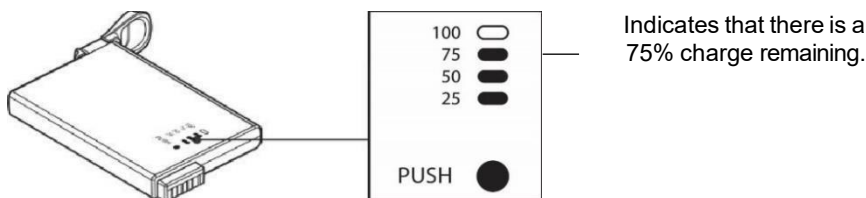


The typical time to recharge a fully discharged battery is approximately 1.5 hours if charging a single battery and 2.5 hours if charging two with the device plugged into the AC power supply and powered off. Battery charge times while the device is running will vary depending on the selected setting.

To maximize battery life:

- Store battery in a cool, dry place when not in use. Store battery with at least 50% capacity remaining.
- If using multiple batteries, uniquely label each battery to ensure each battery is rotated equally.
- Charge batteries every three (3) months when not in use.

There is no routine maintenance or service required for the Oxlife Independence® rechargeable batteries other than the periodic replacement as described above.

AC Power Supply

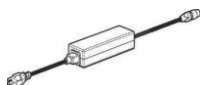

- ❗ Use only the AC power supply (800-1022) provided with this device. Use of a non-approved AC power supply may cause overheating or damage to the device and will void the warranty.
- ❗ The power supply is not water resistant.
- ❗ **DO NOT** place anything in the power supply port other than the supplied AC or DC power cords.
- ❗ Ensure the power supply is in a well-ventilated area. The power supply may become hot during operation. Allow the power supply to cool before handling.

The AC power supply consists of the following components:

- Power supply brick with attached power supply cable to connect to the device.
- AC power input cable

When powered on, the green LED on the power supply will be illuminated.

The AC power supply charges the battery using a 100-240V 50/60 Hz outlet (a typical wall outlet in your home). Using the AC power supply allows you to use your Oxlife Independence[®] while simultaneously charging the battery.

To use the AC power supply, connect the power cord to the AC power converter (brick). Then connect the power supply to a wall outlet and then the Oxlife Independence[®]. The external power icon will be displayed on the control panel.

The AC power supply will charge the battery at all settings.


Recommendation for Use:

Do not use electrical outlets controlled by a switch.

When changing power sources wait for the control panel to display the new power icon and the device to beep before removing the original source.

The detachable power supply cord to the power supply is the means of isolating the device from the supply mains.

If the AC power supply is removed from the device, wait ten (10) seconds before re-applying power.

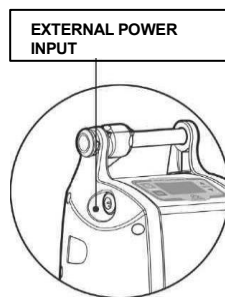

DC Power Cord



The DC power cord allows use of your Oxlife Independence® from a vehicle's 12-volt DC outlet in all settings while the vehicle is running. A 12-volt, 10-amp outlet is recommended for use with the device. DC power is not sufficient to charge batteries at all settings; charging will not occur at settings higher than continuous flow 2.0 or pulse mode 4.0. The external power source icon will be displayed on the control panel.

Each vehicle's DC power outlet varies in specifications and performance. Moreover, such performance may be adversely impacted by one or more applications simultaneously drawing DC power (e.g., cell phones, media players, other electrical systems in the vehicle). In addition, DC power may fluctuate more significantly in stop and go driving conditions.

If a vehicle's DC voltage drops to 11.6 volts or below, the device will alert and revert to battery power. In this condition, the device performance should be monitored; reducing the flow setting may improve the performance of the DC power outlet. If a vehicle's DC voltage drops below 10.6 volts, the device will alarm, and if battery power is not sufficient, the device will shut down. Switch to battery power to continue using your device.

- ❗ Use only the DC power cord (800-1001) provided with this device. Use of power cords not supplied by O2 Concepts® may cause overheating or damage to the device and will void the warranty.

Recommendations for Use:

When operating the Oxlife Independence® in your vehicle, ensure that the device is securely stowed and will not get damaged during transport.

Ensure that the air inlet and exhaust ports are not blocked.

Batteries will not charge at any setting if the engine is not running.

Use the DC outlet closest to the car battery.

The device may "rev" more while operating on DC power. This is normal operation while connected to DC power and is not cause for concern.

DO NOT leave the device plugged into the vehicle when the engine is not running.

DO NOT use the Oxlife Independence® with any power splitting devices.

DC Power Cord User Information Guide

Starting the device on DC power

1. **ALWAYS** have the vehicle's engine running **BEFORE** plugging in your Oxlife Independence®.
2. Plug DC power cord into the device **BEFORE** plugging into the vehicle's DC outlet.
3. Plug DC power cord into the vehicle's DC outlet.
4. Remove the DC power cord from the vehicle's DC outlet when the engine is not running.

Disconnecting the device from DC power

1. **WITHOUT** batteries installed, keep the vehicle's engine running and power the device Off by pressing and holding the power button.
2. When the device is OFF, disconnect the DC power cord from the vehicle's DC outlet.
3. **WITH** batteries installed, disconnect the DC power cord from device. The device will beep once and switch to internal battery power.

DC Power Cord Troubleshooting Guide


Event	Solution
Batteries Not Charging	• Reduce the device setting to the lowest setting on either Pulse or Continuous mode and wait for charging icon to appear. • Slowly increase settings to determine if charging can be achieved at the prescribed setting. Each vehicles' DC power performance varies and may limit achieving the desired flow rate. • Batteries will not charge on settings higher than 2.0 Continuous or 4.0 Pulse
Depleted Battery Charging	• A fully depleted battery may not begin charging unless the device is powered off. The lightning bolt may not display initially in the battery icon until the battery has charged to a minimum level. To prevent accidental discharge of fully charged batteries, batteries can be removed from the device while running on DC power.
Device Alerting Beeping Intermittently	• Occurs when the DC voltage supplied by the vehicle drops in and out of the minimum DC input threshold required for device use. • Using a lower flow setting* may improve DC outlet performance.
Low External Power Alarm	• To reduce potential device shut down on DC power, eliminating other sources of power draw can improve performance. Air conditioning is an example of an electrical system that draws power away from the DC outlet. • Using a lower flow setting* may improve DC outlet performance. • Check your vehicle's user manual or consult an auto technician to determine if your vehicle's DC power amperage and wattage are sufficient. The device requires 150 watts to run at all settings. If available, using a 12V, 15A outlet is preferred. • You may require an inverter with at least a 300- watt capacity. This can be purchased through any auto service provider.

*Consult your healthcare provider for advice on changing settings.

Handle Operation

OPEN/CLOSE

Press button to unlock the handle. Begin rotating and release button. Handle will re-lock in open/close position automatically.

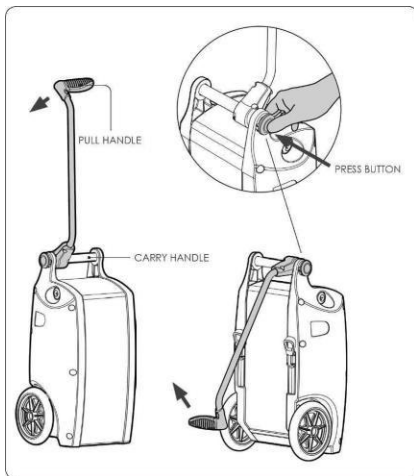

HANDLING

DO NOT LEAN ON THE DEVICE.

To prevent damage to the device, avoid resting excessive weight on the handle.


STAIRS

When transporting the device downstairs, allow the device to travel down first.



Accessory Bag

The Accessory Bag is designed to carry your Oxlife Independence® accessories including batteries, AC power supply with cord, and DC power cord.

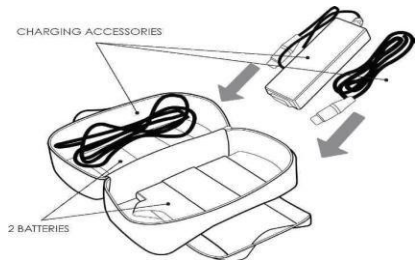

Cannula Use

- ⚠ Your healthcare provider or authorized dealer will recommend the proper cannula for your use.
- ⚠ Ensure selected cannula is compliant to ISO 80601-2-69 Medical Electrical Equipment, particular requirement for basic safety and essential performance of oxygen concentrator equipment.
- ⚠ Use of a cannula not specified for use with oxygen may impair the performance of your device.
- ⚠ **DO NOT** use cannula tubing longer than 7 feet (2m) in Pulse Flow Mode.
- ⚠ **DO NOT** use cannula tubing longer than 50 feet (15m) in Continuous Flow Mode.
- ⚠ A single lumen nasal cannula rated at 6 liters per minute is required to ensure proper oxygen delivery.
- ⚠ To ensure proper oxygen flow, confirm that the cannula is not pinched or blocked in any way.
- ⚠ The cannula fitting may be tight.
- ⚠ **DO NOT** use grease or oils to lubricate the oxygen outlet port.



Read and follow the instructions, included with the cannula and provided by your authorized dealer.

Regularly clean and replace your cannula as instructed by the manufacturer and authorized dealer or healthcare provider.

O2 Concepts recommends the following cannulas by Salter Labs:

Pulse Mode Use: 7ft (2m) cannula: Part #16SOFT-7, or equivalent.

Continuous Mode Use: 50ft (15m) cannula: Part #OTC-203-50, or equivalent.

Firesafe™ Cannula Valve

The Firesafe™ Cannula Valve is provided in order to comply with IEC 60601-1 and ISO 80601-2-69 requirements.

Intended Use:


The Firesafe™ Cannula Valve is a thermal fuse designed to extinguish an oxygen delivery tube fire and stop the flow of oxygen if the tube is accidentally ignited.

- ⚠ Read Firesafe™ Cannula Valve Instructions for Use for the RES010KIT from Sunset Healthcare Solutions before installing. As with all medical equipment, attempting to use or install this device without a thorough understanding of its operations and limitations may result in patient or user injury.

Accessing the Provider Screen


The Provider Screen shows information about the status and language setting of the device. Within the Provider Screen, the language of the device can be changed to six different languages, including English, German, Dutch, Spanish, French, and Italian.

Entering the Provider Screen:

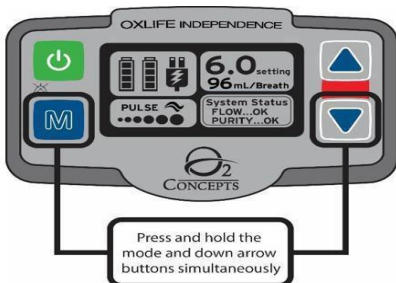

Press and hold the Mode Button  and Down Arrow Button simultaneously to enter the Provider Screen.

Changing Languages

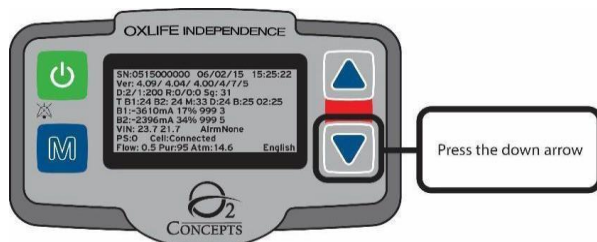

Device Language can be changed by the user.

To change device language:


From the Provider Screen, the user can change the device language by pressing the Down Arrow Button . The current language is displayed in the bottom right corner. The language changes each time the Down Arrow Button  is pressed.

Cell Status and Signal Strength:



The cell signal strength is listed in the third row on a scale of 0 through 31 (ex. Sg:31). The cell status is listed in the eighth row (ex. Cell: Connected).

Voltage In:

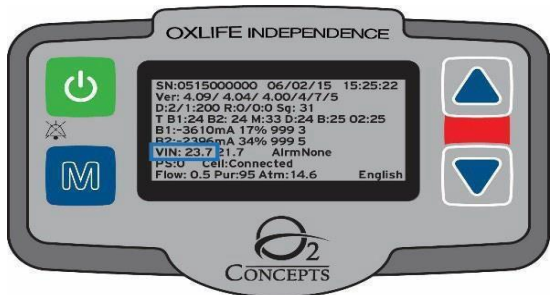


The voltage in is the voltage that the device is receiving from an external power source, such as the AC power supply or DC power cord (ex. VIN: 23.7 where 23.7 is in volts).

Section 5: User Alerts and Alarms

The functionality of the alarm system is verified automatically by the device upon start up by the display flashing the three (3) visual display colors and audible indicators. There will be delays that are greater than ten (10) seconds inherent to specific alarms, and these delays are explained in detail for each alarm in the following tables. To ensure that audible alarm notifications are heard, the maximum distance the user can be from the device is dependent on surrounding noise levels. When two alarm conditions occur at the same time, the higher priority condition will be displayed.

Alarm and Alert Screens

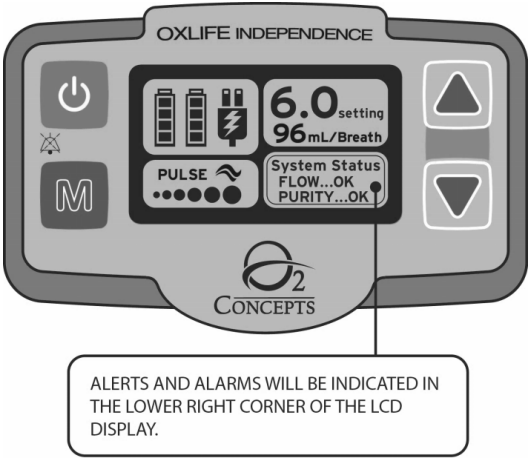

Indefinite Acknowledge (Alarm Mute)

In an alert or alarm condition, pressing the power button momentarily will put the device into an indefinite acknowledge state or alarm mute. In this state, the alarm will be silenced, and the screen will stop flashing. The device will stay in this state until powered off or until a higher-priority alarm is activated.

DO NOT IGNORE ALERTS OR ALARMS

Alarm and Alert Screen Descriptions

Alarms

High Priority alarms are indicated by a red screen flashing at 2Hz, and a 10-beep sequence repeated after 2.5 seconds of silence. All alarms are technical and high priority.

The device will continue to sound the alarm until the power button  is pressed momentarily to silence the audio; the power button is held to power off the device or 15 seconds has expired at which time the device will power off.

Low Purity Alarm	When oxygen levels have drop below 72% for 60 seconds, the text "Low Purity" will flash in the system status window.
Solutions:	Clean or replace air inlet filter. Change to another source of oxygen and contact your authorized dealer.
No Flow Alarm	When the flow of oxygen has stopped for 45 seconds, the text "No Flow" will flash in the system window for up to 30 minutes. If at the end of 30 minutes the No Flow condition has not been corrected, the device will initiate a No Flow Alarm and shut down after 15 seconds. If the No Flow condition is corrected, the device will return to normal function.
Solutions:	Check the cannula connection. Repair or replace tubing. Clean or replace air inlet filter. Move device to ensure adequate air flow. Change to another source of oxygen and contact your authorized dealer.
Low Battery Alarm	When battery power is depleted, the pump shuts off and the text "Low Battery" will flash in the system status window.
Solutions:	Plug into an external power source. Replace depleted battery (batteries) with a charged battery.
High External Power Alarm	When the voltage from an external power source is above 26 volts (as measured by the device internally), the text "High External Power" will flash in the system status window.
Solutions:	Connect to an authorized O2 Concepts Oxlife Independence power supply.
Low External Power Alarm	When the voltage from an external power source falls below 10.6 volts (as measured by the device internally, the text "Low External Power" will flash in the system status window.
Solutions:	Ensure all charging connections are intact. Ensure the battery is installed. Refer to the DC user guide section of this manual. Change to another source of power and contact your authorized dealer.
Unauthorized Battery Alarm	If the device detects that a battery other than one from O2 Concepts has been installed in the device, the battery icon will show an exclamation point within the battery outline. If no other power source is present, the text "Unauthorized Battery" will flash in the system status window.
Solutions:	Remove the unauthorized battery and replace it with an authorized O2 Concepts Oxlife Independence battery.

Over Temperature Alarm	The maximum operating temperature of the enclosure and pump is 70°C and 90°C respectively. If either of these component temperatures are reached for 15 seconds, the text "Over Temperature" will flash in the system status window.
Solutions:	Relocate the device to improve airflow or to a cooler environment and allow the device to reach operating temperature. Utilize an alternate source of oxygen if necessary and contact your authorized dealer.
Invalid Motor Temperature Alarm	When the motor temperature is out of a valid range for 15 seconds the text "Invalid Motor Temperature" will flash in the system status window.
Solutions:	Relocate the device to improve airflow or to a cooler environment and allow the device to reach operating temperature. Utilize an alternate source of oxygen if necessary and contact your authorized dealer.
Invalid Box Temperature Alarm	When the internal box temperature is out of a valid range for 15 seconds, the text "Invalid Box Temperature" will flash in the system status window.
Solutions:	Relocate the device to improve airflow or to a cooler environment and allow the device to reach operating temperature. Utilize an alternate source of oxygen if necessary and contact your authorized dealer.
No External Power Alarm	If power is removed while the device is running it will sound a 10-beep sequence after 4 seconds of silence at which time the device will power off. Since the power is removed, the display is blank, and no backlight color is available.
Solutions:	Ensure all charging connections are intact. In the case of an extended outage, utilize an alternate, non-powered source of oxygen.

DO NOT IGNORE ALERTS OR ALARMS

Alerts

Medium Priority alerts are indicated by an amber screen and a 3-beep sequence repeated after 15 seconds of silence. The No Flow, Unauthorized Battery, and Low External Power alerts are technical alerts and classified as medium priority. The No Breath alert is a physiological alert classified as medium priority.

Low Priority alerts are indicated by an amber screen and a single beep every 30 seconds. Pushing the power button  once silences the alert. The Low Purity and Low Battery alerts are technical alerts classified as low priority.

Low Purity Alert	When oxygen levels drop below 85% for 60 seconds, the text "Low Purity" will flash in the system status window.
Solutions:	Clean or replace air inlet filter. Contact your authorized dealer.
No Flow Alert	When the flow of oxygen has stopped for 45 seconds, the text "No Flow" will flash in the system window for up to 30 minutes. If the No Flow condition is corrected, the device will return to normal function. If after 30 minutes the "No Flow" Alert condition has not been corrected, the device will initiate a No Flow Alarm.
Solutions:	Check the cannula connection. Repair or replace tubing. Clean or replace air inlet filter. Move device to ensure adequate air flow. Change to another source of oxygen and contact your authorized dealer.
Low Battery Alert	When calculated battery run time reaches 6 minutes remaining, the text "Low Battery" will flash in the system status window.
Solutions:	Connect to an external power source. Replace depleted battery (batteries) with a charged battery.
Unauthorized Battery Alert	If the device detects that a battery other than one from O2 Concepts has been installed in the device, the battery icon will show an exclamation point within the battery outline and the text "Unauthorized Battery" will flash in the system status window. The device will continue to run if on AC, DC, or a second valid battery is present.
Solutions:	Remove the unauthorized battery and replace it with an authorized O2 Concepts Oxlife Independence battery.
Low External Power Alert	When the voltage from an external power source falls to 11.6 volts or less (as measured by the device internally) for 3 seconds, the text "Low External Power" will flash in the system status window.
Solutions:	Ensure all charging connections are intact. Ensure a battery is installed. Refer to the DC user guide section of this manual. Change to another source of power and contact your authorized dealer.

No Breath Alert	When in Pulse mode, if a breath is not detected for 45 seconds, the text "No Breath" will flash in the system status window. If no breath is detected for an additional 15 seconds, the device will switch to the previous CONTINUOUS mode setting selected. Note: No Breath Alert is used when in Pulse Mode ONLY.
Solutions:	Check the cannula for kinks or obstructions. Ensure the patient is breathing through their nose. Ensure cannula tubing does not exceed 7 feet (2.1m).


DO NOT IGNORE ALERTS OR ALARMS

Alarm System Test

Use this test to verify that the alarm system is working properly. This test includes screen color, screen flashing sequence, and audible alarm sequence. Follow the steps in the procedure listed below.


Step 1:

Press and hold both the Up-Arrow button and Down Arrow button simultaneously to enter the Service Mode screen.

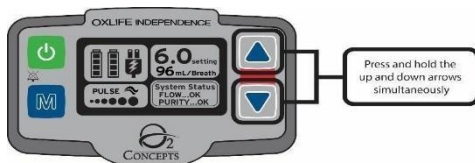

Step 2:

Press the Mode button to test the alarm system.

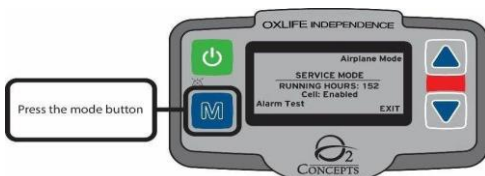

The system will begin by testing the High Priority Alarm. The High Priority Alarm consists of a 10-beep sequence separated by 2.5 seconds of silence, and the LCD screen flashes red.

The device will automatically enter the Medium Priority Alarm test once the High Priority Alarm sequence has been concluded. Pressing the Mode button while the High Priority Alarm test is in progress will terminate the High Priority Alarm test and start the Medium Priority Alarm. The Medium Priority Alarm consists of a 3-beep sequence separated by 15 seconds of silence, and the LCD screen will flash amber.

The device will automatically enter the Low Priority Alarm test once the Medium Priority Alarm sequence has concluded. Pressing the Mode button while the Medium Priority Alarm test is in progress will terminate the Medium Priority Alarm test and start the Low Priority Alarm. The Low Priority Alarm consists of 1 beep and a solid amber LCD screen.

At the conclusion of the test, the screen will display Test Done. Press the Mode button once to reset and complete the Alarm Test feature.

Press the Down Arrow button to exit the Service Mode screen.

Section 6: Maintenance and Cleaning

Maintenance

Routine cleaning of the air inlet filter, and the device care and cleaning, as described below, are the only routine maintenance by the user necessary for the operation of the Oxlife Independence®.

There is no routine maintenance or service to the Oxlife Independence® lithium-ion rechargeable batteries other than periodic replacement as described in the section on Battery Time Management, page [24](#).

All other maintenance or service **MUST** be conducted by a qualified Oxlife Independence® service technician. **DO NOT** attempt to disassemble or perform any maintenance on your device. Any such attempt will void the warranty.

Pre-Use Functional Check

A trained service technician can qualitatively verify that the device is ready for use by performing the following protocols:

- Connect the nasal cannula to the cannula port on the device. Set the device to Continuous flow mode. The technician should be able to hear the flow of oxygen by placing the outlet of the nasal cannula near his or her ear (do not direct flow straight at ear to avoid discomfort).
- The technician should be able to feel the flow of oxygen by placing a finger roughly a half inch from the outlet of the nasal cannula.
- Place the end of the nasal cannula in a half-full cup of water. Verify that there are bubbles from the output of the nasal cannula.


Device Care and Cleaning

Air Inlet Filter Cleaning

Routine cleaning of the air inlet filter as described below and periodic cleaning of the exterior of the device is the only routine maintenance.

To keep your Oxlife Independence[®] working properly, it is recommended you clean the air inlet filter weekly. If the Oxlife Independence[®] is used in a dusty environment, more frequent cleaning of the air inlet filter may be required.

Follow these simple directions to clean the air inlet filter:

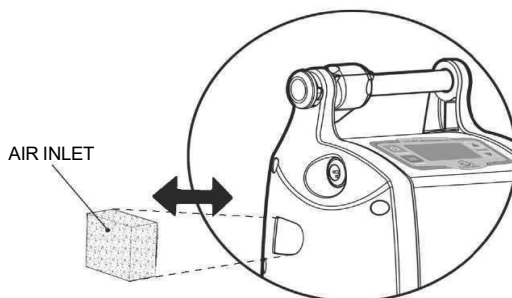

1. Remove the air inlet filter.
2. Wash the filter by running under warm tap water using a mild detergent.
3. Rinse thoroughly under warm, running water.
4. Squeeze out excess water.
5. It is important to allow the filter to dry completely before reinserting it into the device. Moisture from a wet or damp air inlet filter may damage your device.
 - ⚠ Never use your Oxlife Independence[®] without an air inlet filter.
 - ⚠ We recommend that you keep an extra air inlet filter to use as a replacement while one is drying. Contact your authorized dealer for extra air inlet filters.

Device and Battery Cleaning

- ⚠ Unplug your Oxlife Independence before cleaning.
- ⚠ Do not submerge your device or any accessories or allow water to enter the device.

Clean the exterior of your device using a soft cloth dampened with a mild detergent and water and wipe dry. Clean the outside of the device monthly or more frequently, as needed.

- ⚠ **DO NOT** spray or soak the case or front panel.
- ⚠ **DO NOT** spray or soak battery contacts.
- ⚠ **DO NOT** use alcohol, ethylene chloride or petroleum- based cleaners on the case or power supplies.

AC Power Supply and DC Cord Cleaning

- ❗ Unplug your Oxlife Independence[®] and/or power supply before cleaning.

Clean the power supplies and cords using a soft cloth dampened with a mild detergent and water and wipe dry.

Cannula Cleaning

Clean and replace the cannula as instructed by the manufacturer and your healthcare provider.

Disinfection

Between patients, the exterior of the device and battery should be disinfected.

Clean the device as instructed above.

Per manufacturer's instructions, disinfect exterior and battery using Oxycide (or similar product) with a soft cloth.

- ❗ DO NOT use alcohol, ethylene chloride, benzalkonium chloride or petroleum-based cleaners.

Device Disposal

This product may contain substances that could be harmful to the environment if inappropriately disposed of in landfills. Follow local governing ordinances and recycling plans regarding disposal of the device.

Battery Disposal

The battery is rechargeable and can be recycled. Always return batteries to your healthcare provider or an authorized dealer for proper disposal. You can also contact your local government offices for proper disposal instructions for a lithium-ion battery.

Section 7: Traveling

Before traveling, be sure to pack the following:

- › Accessory Bag with User Manual
- › AC Power Supply and DC Power Cord
- › Fully Charged Battery (and extra batteries if required)

Also, bring contact information for your healthcare provider or authorized dealer.

When traveling internationally, a standard, applicable international power plug adapter is required.

Please note that many air carriers require prior notification (usually at least 48 hours) for passengers traveling with portable oxygen concentrators. In addition, most carriers will also require advanced submission of a medical verification statement from your **healthcare provider**. You should confirm any notification and verification requirement with your specific carrier.

Not all air carriers provide an electrical outlet aboard the aircraft, so you should have sufficient batteries (or confirmed alternative power source) to account for at least 150% of your total travel time, including, but not limited to: (A) commuting time to and from the airport; (B) transition time in and out of airports; (C) the duration of your flight; and (D) any unexpected delays. Please be aware that your battery requirements may vary based upon your personal setting.

If traveling by train, bus or boat contact your carrier to inquire about power port/outlet availability.

Entering and Exiting Airplane Mode

The cellular connection must be deactivated prior to flight on any commercial aircraft.

Step 1:

Press and hold both the Up-Arrow button and Down Arrow button simultaneously to enter the Service Mode screen.

Step 2:

Press the Up-Arrow button to enter Airplane Mode. The mode listed in the upper right corner of the display is the mode that you will be entering.

Step 3:

The device will exit Airplane Mode after 24 hours have passed or the user manually re-enters Normal Mode.

Press the Up-Arrow button again to enter Normal Mode.

Step 4:

Press the Down Arrow button to exit the service mode screen.

Section 8: Standards Compliance

This device is designed to comply with the following standards:

- IEC 60601-1: 2005 Ed.3+A1; A2 Medical Electrical Equipment; Part 1: General Requirements for Basic Safety and Essential Performance
- IEC 60601-1-2: 2020 Ed. 4.1; Medical Electrical Equipment, Part 1-2: General Requirement for Safety – Collateral Standard: Electromagnetic (EMC) Compatibility
- IEC 60601-1-6:2010 Ed3+A1; A2 Medical Electrical Equipment – Part 1-6 General Requirements for Basic Safety and Essential Performance – Collateral Standard: Usability
- IEC 60601-1-8:2006 Ed.2 + A1; A2 Medical Electrical Equipment – Part 1-8: General Requirements for Basic Safety and Essential Performance – Collateral Standard: General Requirements Tests and Guidance for Alarm Systems in Medical Electrical Equipment and Medical Electrical Systems
- IEC 60601-1-11:2015 Ed.2+A1 Medical Electrical Equipment Part 1-11: General Requirements for Basic Safety and Essential Performance Collateral Standard – Requirements for Medical Electrical Equipment and Medical Electrical Systems Used in the Home Healthcare Environment
- IEC 62304:2006 Ed1 + A1 Medical Device Software – Software Life Cycle Processes
- ANSI/AAMI IEC 62366-1:2015 Ed.1+A1 Medical Devices – Application of Usability Engineering to Medical Devices
- ISO 80601-2-67:2020 Ed.2 Medical Electrical Equipment – Part 2-67 Particular Requirements for Basic Safety and Essential Performance of Oxygen Conserving Equipment
- ISO 80601-2-69:2020 Ed.2 Medical Electrical Equipment – Part 2-69 Particular Requirements for Basic Safety and Essential Performance of Oxygen Concentrator Equipment
- RTCA, DO 160, Section 21, Category M; Emission of Radio Frequency Energy
- ISO 13485:2016 Medical Devices; Quality Management Systems; Requirements for Regulatory Purposes
- ISO 15223-1:2021/Amd.1:2025 Medical Devices – Symbols to be used with medical device labels, labelling and information to be supplied by the manufacturer – Part 1: General Requirements

Classification

- › IEC Class II Internally Powered Equipment
- › Type BF Applied Part
- › IP22: Drip Proof Equipment (Protected against solid objects over 12 mm and direct sprays of water up to 15° of vertical per IEC 60529)
- › **NOT** suitable for use in the presence of flammable anesthetic mixture with air or with oxygen or nitrous oxide
- › Continuous Operation

OXLIFE INDEPENDENCE®



O2 Concepts, LLC
199 Park Road Extension
Middlebury, CT 06762 USA
+1(877) 867-4008
www.o2-concepts.com



MED IVD HEALTHTECH SL
C/Alejandro Casona, 2 portal 6-3 b 2,
28702 San Sebastian de los Reyes,
Madrid, Spain
+34 93 609 62 69
Email: ecrep@mdxcro.com



2797



© 2026 O2 Concepts, LLC, all rights reserved; O2 Concepts®, Oxlife Independence®, and DNA Technology®, including certain related logos, are registered trademarks of O2 Concepts, LLC and may not be replicated without express written consent.

PN: 800-1097 Rev C – 2026-05-27